

Product Data Sheet

Upright 16 cu. ft. Solid Door Refrigerator, High Performance - Certified to NSF/ANSI 456 Standard for Vaccine Storage

PH-ABT-NSF-16S

Product Description

These cutting-edge pharmacy refrigerators are certified in accordance with the NSF/ANSI 456 Standard for Vaccine Storage. With this certification, units protect pharmaceuticals at optimal temperatures, preventing waste and allowing for peak delivery. Our Premier line includes premium features such as extensive alarm systems and digital touch pad displays.

These solid door refrigerators utilize microprocessor controllers and feature temperature alarms, remote alarm contacts, and probe access ports with included probes. Units run on natural, hydrocarbon refrigerant for environmental health and energy efficiency.

General Description and Application			
Description	Single Solid Door Pharmacy/Vaccine Upright Refrigerator		
Operational environment	Indoor use only, +18°C to +26°C (+65°F to +78°F), <70% RH		
Storage capacity	16 cu. ft. gross volume		
Door	One swing solid door, self-closing, right hinged, non-reversible, magnetic sealed gasket, keyed lock		
Shelves	Six shelves (five adjustable/one fixed) with guard rail on back		
Mounting	3 1/2" Swivel Casters(two locking)		
Interior lighting	Shielded, switched LED lighting, full coverage, balanced spectrum		
Airflow management	Forced Air technology, patent pending		
External probe access	Rear wall port (3/4") dia.		
Insulation	Cabinet is foamed-in-place with EPA compliant high density urethane foam		
Exterior materials	White powder coated steel		
Access control	Pyxis [®] , Omnicell [®] and AcuDose RX [®] compatible		
General warranty	Two (2) years parts and labor warranty, excluding display probe calibration		
Compressor warranty	Five (5) years compressor warranty		
Product Weight	243 lbs.		
Shipping Weight	283 lbs.		
Rated Amperage	3 Amps		
Power Plug/Power Cord	NEMA 5-15 plug, 8 to 10 ft typical, conforms to UL471 requirements, Vaccine Storage power cord warning label		
Facility Electrical Requirement	110-120V AC: 15 A (minimum)		
Agency Listing and Certification	Certified with the temperature performance requirements as defined in the NSF/ANSI 456 Standard for Vaccine Storage for all testing scenarios. UL, C-UL, ETL, C-ETL listed and certified to UL471 standard, hydrocarbon refrigerant safety. Energy Star Certified		
Included Accessories	Temperature monitor device (TMD) complies with the current CDC guidelines, with 3 years certification of calibration, "buffered" probe in the product simulated solution, min/max memory, field installable, and visual & audible temp alarm		
	Pharmacy refrigerator/freezer toolkit and temperature logs		


Refrigeration System	
Compressor	Hermetic, high performance
Refrigerant	EPA SNAP compliant, R290, propane
Condenser	Fin and tube design, high efficiency fan
Evaporator	Fin and tube design, high efficiency fan
Defrost	Cycle optimized, zero energy

Certifications




*-one or more of these certifications may apply to this unit.

Temperature Probes					
Probe	Ave	Min	Max		
1	4.3	3.4	5.4		
2	4.6	4.2	5.2		
3	4.7	4.4	5.2		
4	4.6	4.1	5.4		
5	4.6	4.2	5.3		
6	4.3	3.7	5.3		
7	4.6	4.3	5.2		
8	4.7	4.3	5.4		
9	4.7	4.1	5.6		
10	4.6	4.2	5.3		
11	4.6	3.8	5.7		
12	4.5	3.9	5.3		
13	4.7	4.4	5.1		
14	5.4	5.1	5.9		
15	4.3	3.5	5.4		

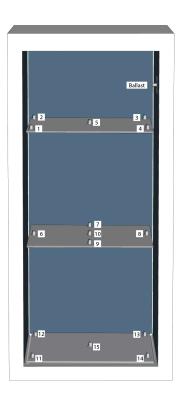

Temperature Charts

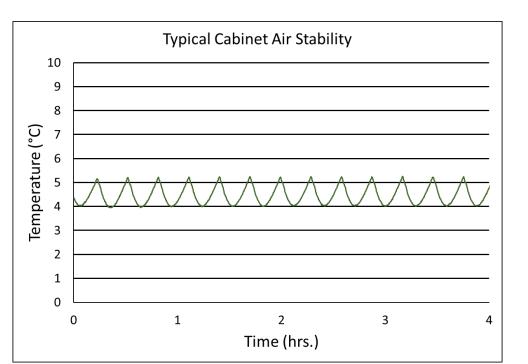

Performance	
Uniformity ¹ (Cabinet air)	+/- 0.9°C
Stability ² (Cabinet air)	+/- 1.0°C
Maximum temperature variation (Cabinet air)	+/-1.2°C
Temperature rise after 8 sec door openings	Temperature did not exceed 6.5°C at any probe for all required NSF/ANSI 456 testing
	protocols ³
Recovery after 3 min door opening	All probes recover to under 8°C within 4.3 min.
Energy consumption	1.25 KWh/day ^₄
Average heat rejection	1.97 KWh/day (280BTU/h)⁴
Noise pressure level (dBA)	48 or less installed
Pull down time to 4°C nominal operating temp	30 min

Controller, Configuration, Alarms and Monitoring

Controller technology	Parametric, microprocessor, LED display with 0.1°C resolution
Display technology	NSF/ANSI 456 Standard for Vaccine Storage compliant digital temperature display and alarm module with battery back-up, F/C switchable.
Temperature setpoint range	1°C to 10°C (Controller settings must remain unaltered to ensure thermal performance compliant with NSF/ANSI 456 Standard for Vaccine Storage requirements)
Display probe	Calibrated, stainless steel
External alarm connection	State switching remote alarm contacts
	Visual and audible indicators
Alarms	High / Low temperature, compliant with alarm requirements defined in the NSF/ANSI 456
	Standard for Vaccine Storage
Simulator ballast	Glass bead thermal media


Performance data acquired at 22°C ambient, using NSF/ANSI 456 compliant validation ballast probes, empty chamber, during stabilized steady state operation and a DAQ sampling rate of one measurement every 10 seconds


1 - Uniformity is defined as the maximum variance in temperature across all probes at any point in time over the testing period

2 - Stability is defined as the maximum variance in temperature experienced by any single probe over the testing period

3 - Temperature performance for all loaded and unloaded door opening protocols, all alarm, controller and probe requirements as defined in the NSF/ANSI 456 standard for vaccine storage

4 - Data per Energy Star test results or equivalent testing and calculation. Heat rejection based on daily averages, not continuous operation. Performance exceeds Energy Star requirements.

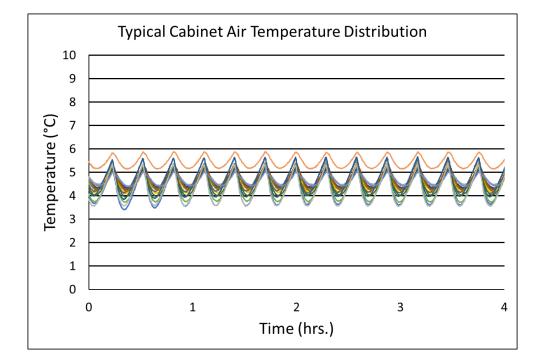

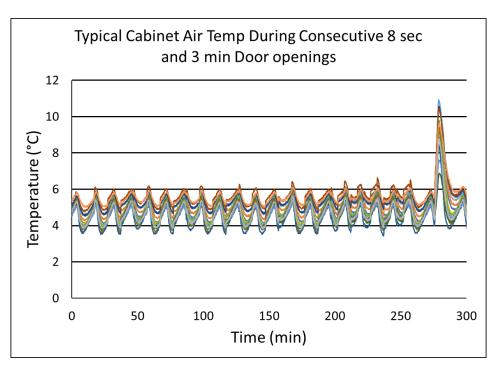



Product Data Sheet Upright 16 cu. ft. Solid Door Refrigerator, High Performance - Certified to NSF/ANSI 456 Standard for Vaccine Storage

Images



Dimensio	Dimensions						
		Width	Depth	Height	Door Swing	Total open Depth	
	Exterior	25"	26"	79"	23 1/8"	49 1/4"	
	Interior	20"	19 1/8"	48 1/8"			

↓ 25" →

