

PH-ABT-NSF-UCFS-0504

Product Description

These premier undercounter refrigerators are designed in accordance with the NSF/ANSI 456 Standard for Vaccine Storage. Units protect pharmaceuticals at optimal temperatures, preventing waste and allowing for peak delivery.

These solid door freestanding refrigerators utilize microprocessor controllers and feature temperature alarms, remote alarm contacts, LED interior lighting, and probe access ports with included probes. American Biotech Supply Vaccine Storage Refrigerators utilize HFC-free refrigerant for environmental health and energy efficiency.

General Description and Application

Single Solid Door Pharmacy/Vaccine Undercounter Refrigerator Freestanding Description Indoor use only. Optimal operating range: +18°C to +26°C (+65°F to +78°F), <70% RH Operational environment

5.2 cu. ft. gross volume Storage capacity

One swing solid door, self-closing, right hinged, non-reversible, magnetic sealed gasket, keyed Door

Shelves Three shelves (two adjustable/one fixed) with guard rail on back

Leveling legs. Note: 4" of clearance on all sides must be maintained for adequate ventilation Mounting and Installation

Shielded, switched LED lighting, full coverage, balanced spectrum Interior lighting

Forced Air technology, patent pending Airflow management

External probe access Rear wall port (3/8") dia.

Cabinet is foamed-in-place with EPA compliant high density urethane foam Insulation

Exterior materials White powder coated steel

Pyxis®, Omnicell® and AcuDose RX® compatible Access control

Two (2) years parts and labor warranty, excluding display probe calibration General warranty

Five (5) years compressor warranty Compressor warranty

96 lbs. **Product Weight** 132 lbs. Shipping Weight 1.3 Amps Rated Amperage

Power Plug/Power Cord NEMA 5-15 plug, 8 to 10 ft typical, conforms to UL471 requirements, Vaccine storage power

cord warning label

110-120V AC: 15 A (minimum) Facility Electrical Requirement

Certified in accordance with the NSF/ANSI 456 Standard for Vaccine Storage. UL, C-UL, ETL, C-Agency Listing and Certification

ETL listed (either single or dual agency listings) and certified to UL471 standard, hydrocarbon

refrigerant safety. Energy Star Certified

Temperature monitor device (TMD) complies with the current CDC guidelines, with 3 years

certification of calibration, "buffered" probe in the product simulated solution, min/max

memory, field installable, and visual & audible temp alarm

Pharmacy refrigerator/freezer toolkit and temperature logs

Refrigeration System

Included Accessories

Compressor Hermetic, high performance EPA SNAP compliant, R600a, Isobutane Refrigerant Condenser Tube and grid construction, fanless Evaporator Plate wall

Cycle optimized, zero energy Defrost

Performance

Uniformity¹ (Cabinet air) +/- 1.4°C Stability² (Cabinet air) +/- 1.3°C Maximum temperature variation +/-1.7°C

Temperature rise after 8 sec door Temperature did not exceed 6.5°C at any probe for all required NSF/ANSI 456 testing

openings scenarios3

Recovery after 3 min door opening All probes recover to under 8°C within 6 min.

1.15 KWh/day⁴ Energy consumption

1.67 KWh/day (237 BTU/h)4 Average heat rejection Noise pressure level (dBA) 41 or less installed

Pull down time to nominal operating 42 min

temp

Simulator ballast

Controller, Configuration, Alarms and Monitoring

Controller technology Parametric, microprocessor, LED display with 0.1°C resolution

Temperature setpoint range 1°C to 10°C (Setpoint must remain unaltered from the factory setting to remain compliant with

NSF/ANSI 456 Standard for Vaccine Storage requirements)

Display probe Calibrated, stainless steel

External alarm connection State switching remote alarm contacts

Visual and audible indicators

High / Low temperature, compliant with alarm requirements defined in the NSF/ANSI 456 Alarms

Standard for Vaccine Storage

Glass bead thermal media

Performance data acquired at 22°C ambient, using NSF/ANSI 456 compliant validation ballast probes, empty chamber, during stabilized steady state operation and a DAQ sampling rate of one measurement every 10 seconds

- 1 Uniformity is defined as the maximum variance in temperature across all probes at any point in time over the testing period
- 2 Stability is defined as the maximum variance in temperature experienced by any single probe over the testing period
- 3 Temperature performance for all loaded and unloaded door opening protocols, all alarm, controller and probe requirements as defined in the NSF/ANSI 456 standard for vaccine storage
- 4 Data per Energy Star test results or equivalent testing and calculation. Heat rejection based on daily averages, not continuous operation. Performance exceeds Energy Star requirements.

Product Data Sheet

Undercounter 5.2 cu. ft. Solid Door Freestanding Vaccine Refrigerator - Certified to NSF/ANSI 456 Standard for Vaccine Storage

Certifications












*-one or more of these certifications may apply to this unit.


| Temperature Probes | | | | | | | |
|--------------------|-----|-----|-----|--|--|--|--|
| Probe | Ave | Min | Max | | | | |
| 1 | 3.9 | 2.7 | 5.2 | | | | |
| 2 | 4.3 | 4.0 | 4.8 | | | | |
| 3 | 4.1 | 3.4 | 4.9 | | | | |
| 4 | 3.9 | 2.7 | 5.2 | | | | |
| 5 | 4.3 | 4.0 | 4.8 | | | | |
| 6 | 4.5 | 4.0 | 5.2 | | | | |
| 7 | 4.8 | 4.4 | 5.3 | | | | |
| 8 | 4.7 | 4.2 | 5.2 | | | | |
| 9 | 4.9 | 4.4 | 5.5 | | | | |
| 10 | 5.1 | 4.6 | 5.6 | | | | |
| 11 | 5.5 | 5.0 | 6.1 | | | | |
| 12 | 5.5 | 5.2 | 5.8 | | | | |
| 13 | 5.5 | 5.3 | 5.8 | | | | |
| 14 | 5.0 | 4.4 | 5.8 | | | | |
| 15 | 5.3 | 5.0 | 5.7 | | | | |

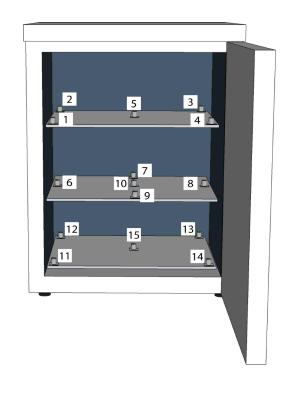

Temperature Charts

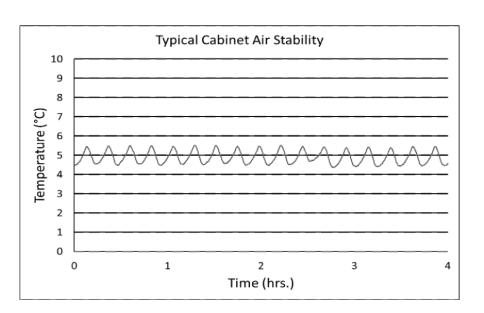

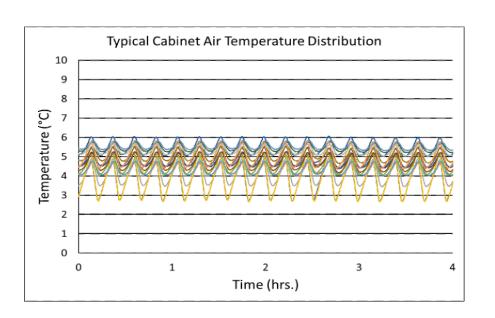

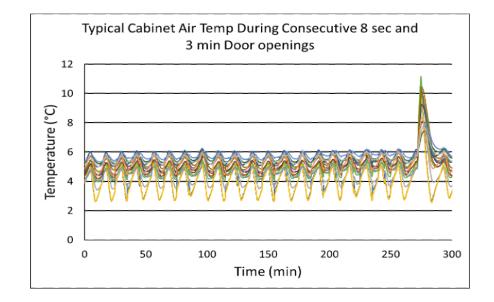



Product Data Sheet


Undercounter 5.2 cu. ft. Solid Door Freestanding Vaccine Refrigerator - Certified to NSF/ANSI 456 Standard for Vaccine Storage

Images





| Dimensions | | | | | | | |
|------------|---------|---------|---------|------------|------------------|--|--|
| | Width | Depth | Height | Door Swing | Total open Depth | | |
| Exterior | 23 3/4" | 24" | 32 1/8" | 22" | 45 3/4" | | |
| Interior | 19 7/8" | 15 3/8" | 24 3/4" | | | | |

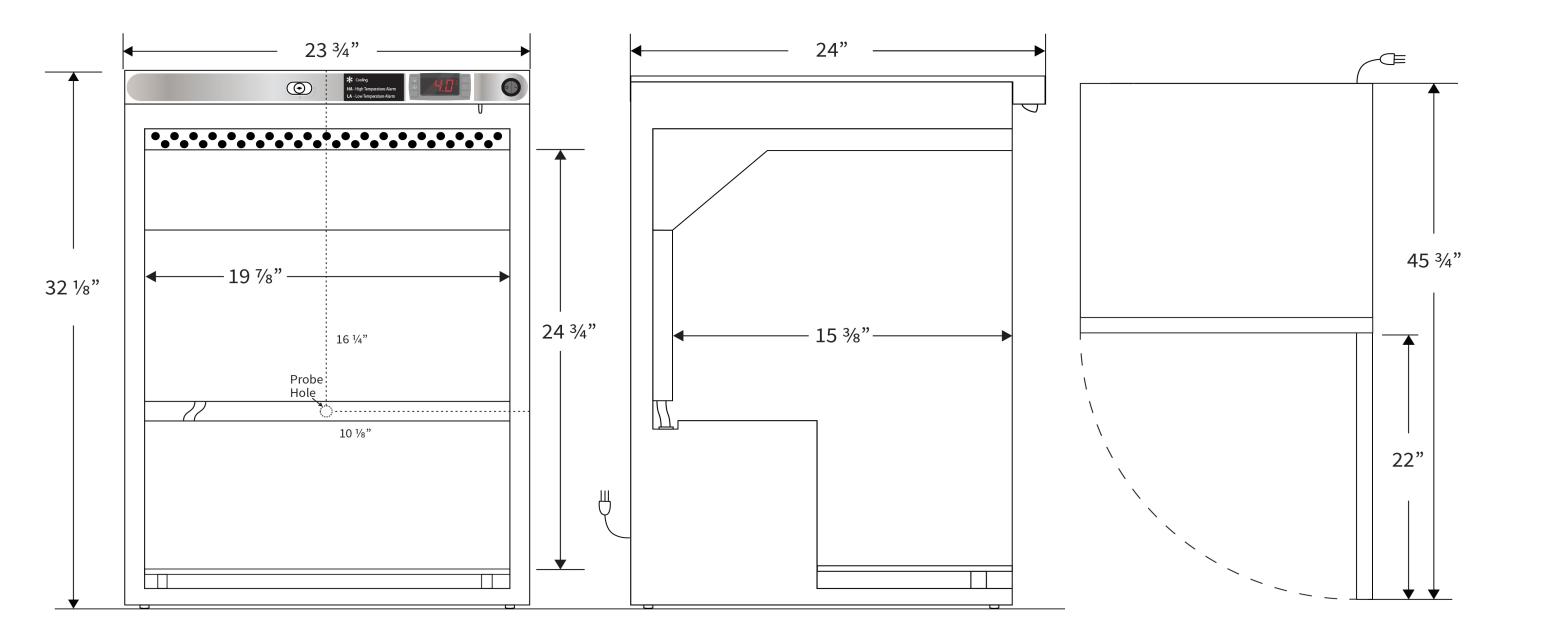

Note: This unit must have 4" clearance on sides and back for adequate ventilation

Rev_10042022