

Product Data Sheet

7 CF Premier Pharmacy/Vaccine Refrigerator & Freezer Combination

PH-ABT-HC-RFC7

Images



Product Description

The Premier Pharmacy Dual Temp Combination Refrigerator/Freezers feature separate microprocessor temperature controller allowing for precise temperature control, verification and recovery. An interior fan driven, forced-draft air circulation system provides excellent temperature uniformity and recovery after door openings REFRIGERATORS ONLY.

All models are equipped with digital temperature displays for convenient monitoring. Critical samples and supplies are security backed with an extensive package of audible and visual alarms. The freezer units feature either an automatic or manual defrost operating system. For manual defrost freezers, periodic manual defrosting is required.

*Please note this unit requires 2 outlets.

Models require 4" of clearance around back and sides of the unit for proper ventilation. Improper installation will void the warranty and will lead to temperature maintenance issues, overheating and possible failure.

Certifications



General Description and Application

Storage capacity (cu. ft)	7 Cu. Ft. capacity (5.2 Ref./1.7 Freezer)
Door	2 exterior swing doors (glass door refrigerator/solid door freezer), right hinged
Shelves	3 adjustable shelves refrigerator/2 adjustable shelves freezer
Freezer Compartments	1
Baskets	Non-applicable
Mounting	Leveling Legs and stacking kit
Interior lighting	Interior light, refrigerator
Airflow management	Forced draft air circulation, refrigerator
External probe access	Two probe access ports (3/4"); one, refrigerator/one, freezer
Insulation	High density urethane foam cabinet insulation
Exterior materials	Steel
Access Control	Keyed door lock
General warranty	Two (2) year parts and labor warranty
Compressor Parts Warranty	Five (5) year compressor parts warranty
Product Weight (lbs)	135
Shipping Weight (lbs)	207
Rated Amperage	1 Amps
Power Plug/Power Cord	NEMA 5-15 Plug, conforms to UL471 requirements . Vaccine storage power cord warning label
Facility Electrical Requirement	110-120V AC: 15 A (minimum)
Agency Listing and Certification	UL-C-UL Listed


Refrigeration System

Refrigerant	Hfc-free refrigerant (r600a)
Defrost	Manual defrost freezer / cycle defrost refrigerator

Controller, Configuration, Alarms and Monitoring

Controller technology	Microprocess temperature controller
Calibration	Calibrated using a NIST traceable device, certificate included
Adjustable Temperature Range	Refrigerator: 36°F – 46°F (2°C -8°C); Freezer: 5°F – -13°F (-15°C - -25°C)
External alarm connection	Remote alarm contacts
Alarms	Audible and visual high and low temperature alarms
Other	Pyxis®, Omnicell® and AcuDose RX® compatible
Data Transfer	Non-applicable
Disclaimers	Models require 4" of clearance around back and sides of the unit for proper ventilation. Improper installation will void the warranty and will lead to temperature maintenance issues, overheating and possible failure. Periodic defrosting of the freezer is required. Product requires 2 outlets.

Dimensions					
	Width (in.)	Depth (in.)	Height (in.)	Door Swing (in.)	Total open Depth (in.)
Exterior	23.625"	24"	53"	22"	46"
Interior	17.75" upper / 19.875" lower	15.75" upper / 16" lower	13.75" upper / 26" lower		

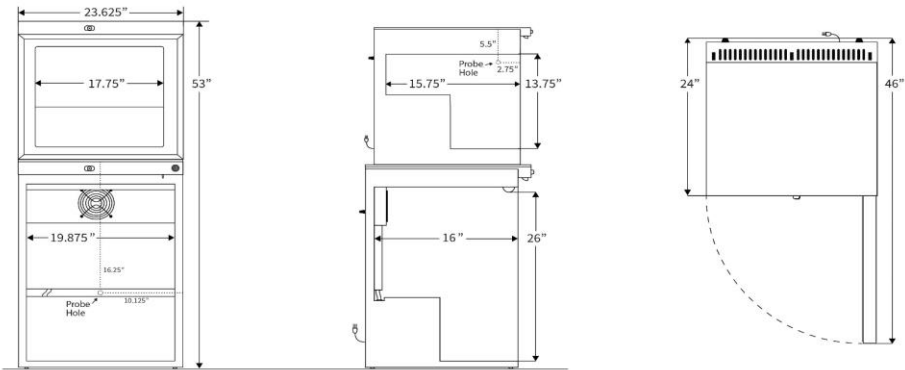

Contact		
Customer Service	800-648-4041 Option 3	customerservice@horizontscientific.com
Technical Service	800-648-4041 Option 5, Parts Option 4	technicalservice@horizontscientific.com