

Product Data Sheet

7 CF Plus Series Pharmacy Refrigerator & Auto Defrost
Freezer Solid Door Combination

PH-ABT-HC-RFC7SA

Images



Product Description

The Plus Series Pharmacy Dual Temp Combination Refrigerator/Freezers feature separate microprocessor temperature controller allowing for precise temperature control, verification and recovery. An interior fan driven, forced-draft air circulation system provides excellent temperature uniformity and recovery after door openings REFRIGERATORS ONLY.

Critical samples and supplies are security backed with an extensive package of audible and visual alarms.

*Please note this unit requires 2 outlets.

Models require 4" of clearance around back and sides of the unit for proper ventilation. Improper installation will void the warranty and will lead to temperature maintenance issues, overheating and possible failure.

Certifications



General Description and Application

Storage capacity (cu. ft.)	5.2 cu. ft. capacity (Refrigerator) / 1.3 cu. ft. capacity (Freezer)
Door	2 exterior swing doors (solid door refrigerator/solid door freezer), right hinged, not reversible
Shelves	3 adjustable shelves refrigerator/2 adjustable shelves freezer
Freezer Compartments	1
Baskets	Non-applicable
Mounting	Leveling Legs and stacking kit
Interior lighting	LED interior light
Airflow management	Forced draft air circulation, refrigerator
External probe access	Two probe access ports (3/4"); one, refrigerator/one, freezer
Insulation	High density urethane foam cabinet insulation
Exterior materials	Steel
Access Control	Keyed door lock
General warranty	Two (2) year parts and labor warranty
Compressor Parts Warranty	Five (5) compressor parts warranty

Product Weight (lbs)	135
Shipping Weight (lbs)	215
Rated Amperage	1 Amps
Power Plug/Power Cord	NEMA 5-15 Plug, conforms to UL471 requirements
Facility Electrical Requirement	110-120V AC: 15 A (minimum)
Agency Listing and Certification	UL-C-UL Listed

Refrigeration System


Refrigerant	HFC-free refrigerant (R600a)
Defrost	Cycle defrost refrigerator and auto defrost freezer

Controller, Configuration, Alarms and Monitoring

Controller technology	Microprocessor temperature controller, refrigerator and freezer
Calibration	Calibrated using a NIST traceable device, certificate included
Adjustable Temperature Range	Refrigerator: 36°F – 46°F (2°C -8°C); Freezer: 5°F – -13°F (-15°C - -25°C)
External alarm connection	Remote alarm contacts
Alarms	Audible and visual high and low temperature alarms

Other	Pyxis®, Omnicell® and AcuDose RX® compatible
Data Transfer	Non-applicable
Disclaimers	Indoor use only, +18C to +26C (+65F to +78F), <70% RH. Specifications are subject to change. Requires 2 outlets. Auto defrost freezers incorporate an electric heater installed on the evaporator that assist in the removal of the frost and ice. Although the heater helps speed up the frost and ice removal, it can raise the temperatures within the freezer. Models require 4" of clearance around back and sides of the unit for proper ventilation. Improper installation will void the warranty and will lead to temperature maintenance issues, overheating and possible failure.

Dimensions					
	Width (in.)	Depth (in.)	Height (in.)	Door Swing (in.)	Total open Depth (in.)
Exterior	23.75"	24"	53"	22"	46"
Interior	17.75" upper / 19.875" lower	9.125" upper / 16" lower	13.75" upper / 26" lower		

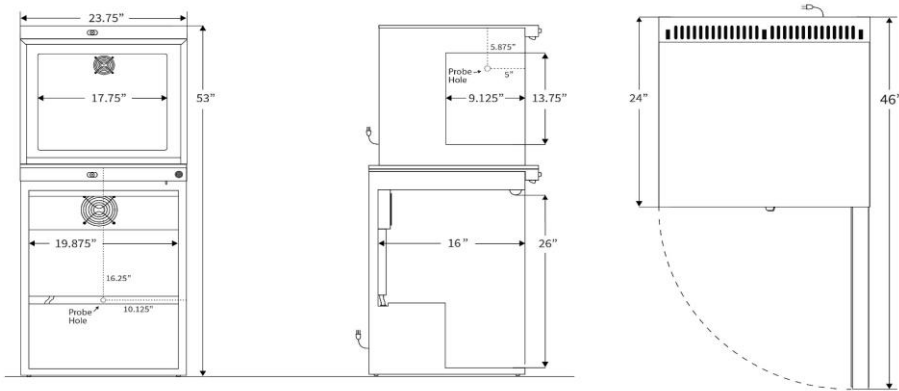

Contact		
Customer Service	800-648-4041 Option 3	customerservice@horizonscientific.com
Technical Service	800-648-4041 Option 5, Parts Option 4	technicalservice@horizonscientific.com