

Product Data Sheet

23 CF Premier Solid Door Stainless Steel Pharmacy Refrigerator

PH-ABT-HC-SSP-23

Images



Product Description

The Premier Pharmacy Stainless Steel Refrigerators feature a microprocessor temperature controller, allowing for precise temperature control, recovery and maintenance versus mechanical style thermostats.

Critical samples and supplies are security backed with an extensive package of temperature security and alarm features.

Environmentally friendly, natural hydrocarbon (HC) refrigerants (not containing hydrofluorocarbons) vastly reduce global-warming potential (GWP) while meeting new EPA/SNAP mandates and UL, ASHRAE compliance guidelines.

Certifications



General Description and Application

Storage capacity (cu. ft)	23 Cu. Ft. Capacity
Door	1 Solid Stainless Steel, Right Hinged
Shelves	3 adjustable
Freezer Compartments	Non-applicable
Baskets	Non-applicable
Mounting	Swivel Casters
Interior lighting	LED Interior light
Airflow management	Forced draft air circulation
External probe access	Probe access port (3/4")
Insulation	High density urethane foam cabinet insulation
Exterior materials	Stainless Steel
Access Control	Keyed door lock
General warranty	Two (2) year parts and labor warranty
Compressor Parts Warranty	Five (5) year compressor parts warranty
Product Weight (lbs)	255
Shipping Weight (lbs)	295
Rated Amperage	6 Amps
Power Plug/Power Cord	NEMA 5-15 Plug, conforms to UL471 requirements. Vaccine storage power cord warning label
Facility Electrical Requirement	110-120V AC: 15 A (minimum)
Agency Listing and Certification	ETL Listed


Refrigeration System

Refrigerant	HFC-free refrigerant (R290)
Defrost	Cycle

Controller, Configuration, Alarms and Monitoring

Controller technology	Microprocessor temperature controller
Calibration	Calibrated using a NIST traceable device, certificate included
Adjustable Temperature Range	36°F – 46°F (2°C- 8°C)
External alarm connection	Remote alarm contacts
Alarms	Audible and visual high and low temperature alarms
Digital Communication	RS-485 (MODBUS)
Warranty Disclaimer	Indoor use only, +18C to +26C (+65F to +78F), <70% RH.
	Specifications are subject to change.

Dimensions					
	Width (in.)	Depth (in.)	Height (in.)	Door Swing (in.)	Total open Depth (in.)
Exterior	26.875"	36.25"	83.25"	24.5"	60.75"
Interior	22"	26.5"	50"		

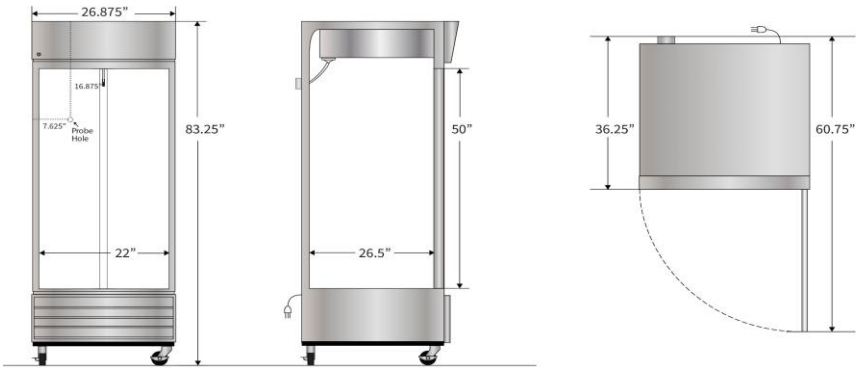

Contact		
Customer Service	800-648-4041 Option 3	customerservice@horizontscientific.com
Technical Service	800-648-4041 Option 5, Parts Option 4	technicalservice@horizontscientific.com

Rev_01202025