

Product Data Sheet

4.6 CF Premier Solid Door Pharmacy/Vaccine
Undercounter Refrigerator Built-In ADA

PH-ABT-HC-UCBI-0404-ADA

Images



Product Description

The Premier Built In Undercounter Refrigerators feature a microprocessor temperature controller and display, which allows for precise temperature control, verification and recovery versus mechanical style thermostats.

Critical samples and supplies are security backed with an extensive package of audible and visual alarms. The Built-In models are designed with a "front breathing" ventilation system, eliminating the need for clearance around the top and sides of the unit; enabling the user to incorporate these units into cabinets and other applications where free space is limited.

Environmentally friendly, natural hydrocarbon (HC) refrigerants (not containing hydrofluorocarbons) vastly reduce global-warming potential (GWP) while meeting new EPA/SNAP mandates and UL, ASHRAE compliance guidelines.

Certifications



General Description and Application

Storage capacity (cu. ft)	4.6 Cu. Ft. Capacity
Door	1 Solid, Right Hinged
Shelves	2 adjustable/1 fixed
Freezer Compartments	Non-applicable
Baskets	Non-applicable
Mounting	Levelling Legs
Interior lighting	Non-applicable
Airflow management	Forced draft air circulation
External probe access	Rear wall port (3/4") dia.
Insulation	High density urethane foam cabinet insulation
Exterior materials	Steel
Access Control	Keyed door lock
General warranty	Two (2) year parts and labor warranty
Compressor Parts Warranty	Five (5) year compressor parts warranty
Product Weight (lbs)	97
Shipping Weight (lbs)	150
Rated Amperage	1.74 Amps
Power Plug/Power Cord	NEMA 5-15 Plug, conforms to UL471 requirements - Vaccine storage power cord warning label
Facility Electrical Requirement	110-120V AC: 15 A (minimum)
Agency Listing and Certification	UL-C-UL Listed, Energy Star Certified, ADA Compliant


Refrigeration System

Refrigerant	HFC-free refrigerant (R600a)
Defrost	Cycle

Controller, Configuration, Alarms and Monitoring

Controller technology	Microprocessor temperature controller
Calibration	Calibrated using a NIST traceable device, certificate included
Adjustable Temperature Range	36°F – 46°F (2°C – 8°C)
External alarm connection	Remote alarm contacts
Alarms	Audible and visual high and low temperature alarms
Other	Pyxis®, Omnicell® and AcuDose RX® compatible
Disclaimers	Indoor use only, +18C to +26C (+65F to +78F), <70%RH Specifications are subject to change.

Dimensions					
	Width (in.)	Depth (in.)	Height (in.)	Door Swing (in.)	Total open Depth (in.)
Exterior	23.75"	24.5"	31.93"	23.5"	46.375"
Interior	19.25"	17.5"	24.375"		

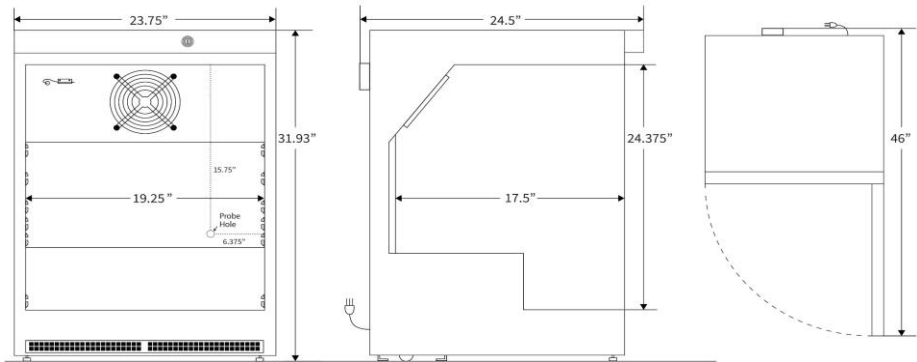

Contact		
Customer Service	800-648-4041 Option 3	customerservice@horizontscientific.com
Technical Service	800-648-4041 Option 5, Parts Option 4	technicalservice@horizontscientific.com

Rev_01202025