

## Product Data Sheet

4.6 CF Plus Series Pharmacy/Vaccine Undercounter Refrigerator Built-In ADA

### PH-ABT-HC-UCBI-0404G-ADA

#### Images



#### Product Description

The Plus Series Built In Undercounter Refrigerators feature a microprocessor temperature controller and display, which allows for precise temperature control, verification and recovery versus mechanical style thermostats.

The Built-In models are designed with a "front breathing" ventilation system, eliminating the need for clearance around the top and sides of the unit; enabling the user to incorporate these units into cabinets and other applications where free space is limited.

Environmentally friendly, natural hydrocarbon (HC) refrigerants (not containing hydrofluorocarbons) vastly reduce global-warming potential (GWP) while meeting new EPA/SNAP mandates and UL, ASHRAE compliance guidelines.

#### Certifications



#### General Description and Application

|                           |                                               |
|---------------------------|-----------------------------------------------|
| Storage capacity (cu. ft) | 4.6 Cu. Ft. Capacity                          |
| Door                      | 1 glass swing door, right hinged              |
| Shelves                   | 2 adjustable/1 fixed                          |
| Freezer Compartments      | Non-applicable                                |
| Baskets                   | Non-applicable                                |
| Mounting                  | Levelling Legs                                |
| Interior lighting         | LED Interior Light                            |
| Airflow management        | Forced draft air circulation                  |
| External probe access     | Probe access port (3/4")                      |
| Insulation                | High density urethane foam cabinet insulation |
| Exterior materials        | Steel                                         |
| Access Control            | Keyed door lock                               |
| General warranty          | Two (2) year parts and labor warranty         |
| Compressor Parts Warranty | Five (5) year compressor parts warranty       |

|                                  |                                                                                           |
|----------------------------------|-------------------------------------------------------------------------------------------|
| Product Weight (lbs)             | 97                                                                                        |
| Shipping Weight (lbs)            | 150                                                                                       |
| Rated Amperage                   | 1.74 Amps                                                                                 |
| Power Plug/Power Cord            | NEMA 5-15 Plug, conforms to UL471 requirements . Vaccine storage power cord warning label |
| Facility Electrical Requirement  | 110-120V AC: 15 A (minimum)                                                               |
| Agency Listing and Certification | UL-C-UL Listed, Energy Star Certified, ADA Compliant                                      |


#### Refrigeration System

|             |                              |
|-------------|------------------------------|
| Refrigerant | HFC-free refrigerant (R600a) |
| Defrost     | Cycle                        |

#### Controller, Configuration, Alarms and Monitoring

|                              |                                                                                               |
|------------------------------|-----------------------------------------------------------------------------------------------|
| Controller technology        | Microprocessor temperature controller                                                         |
| Calibration                  | Calibrated using a NIST traceable device, certificate included                                |
| Adjustable Temperature Range | 36°F – 46°F (2°C – 8°C)                                                                       |
| External alarm connection    | Remote alarm contacts                                                                         |
| Alarms                       | Audible and visual high and low temperature alarms                                            |
| Other                        | Pyxis®, Omnicell® and AcuDose RX® compatible                                                  |
| Disclaimers                  | Indoor use only, +18C to +26C (+65F to +78F), <70%RH<br>Specifications are subject to change. |

| Dimensions |             |             |              |                  |                        |
|------------|-------------|-------------|--------------|------------------|------------------------|
|            | Width (in.) | Depth (in.) | Height (in.) | Door Swing (in.) | Total open Depth (in.) |
| Exterior   | 23.75"      | 25.5"       | 31.93"       | 22.125"          | 46"                    |
| Interior   | 19.25"      | 17.5"       | 24.375"      |                  |                        |

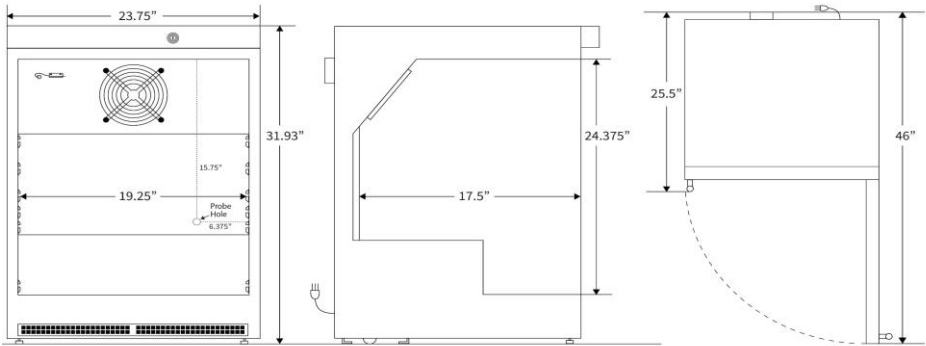

| Contact           |                                       |                                        |
|-------------------|---------------------------------------|----------------------------------------|
| Customer Service  | 800-648-4041 Option 3                 | customerservice@horizonscientific.com  |
| Technical Service | 800-648-4041 Option 5, Parts Option 4 | technicalservice@horizonscientific.com |

Rev\_01202025